



# Live Cell RNA Detection

SmartFlare™ RNA Detection Probes

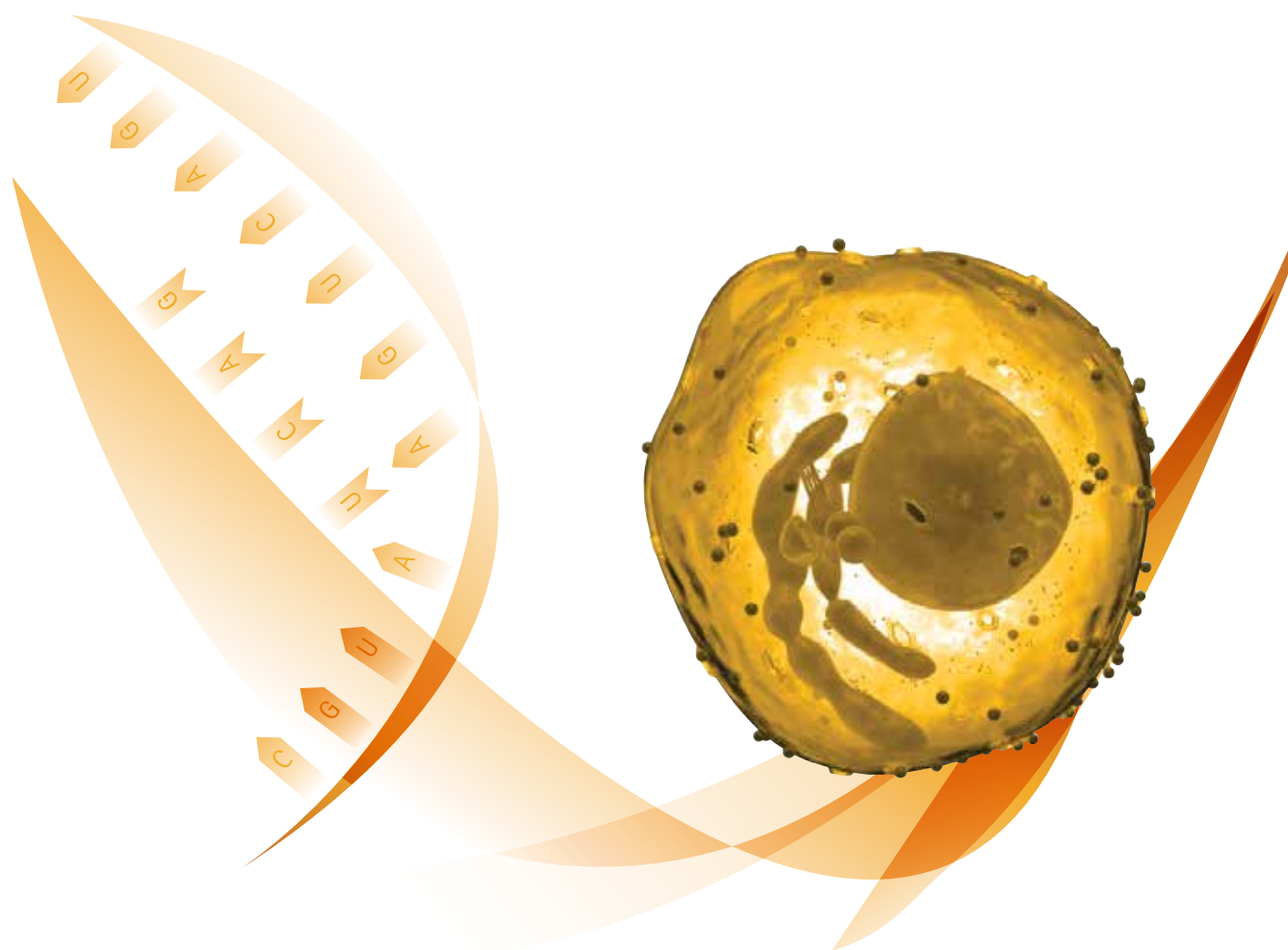

Go live. It's smart.

# Live cell RNA detection with SmartFlare™ RNA Detection Probes

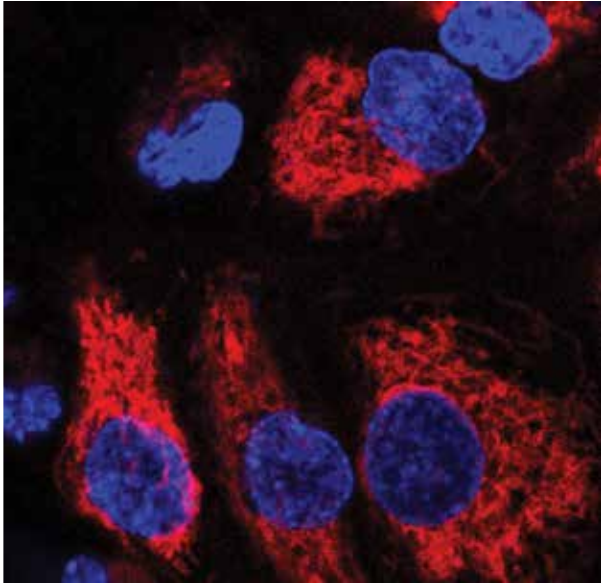

---

## Live cell imaging of miR-21 RNA in human prostate cancer cells (DU145) using SmartFlare™ RNA Detection Probes

DU145 human prostate cancer cells were incubated with 100pM of miR-21 Hu-Cy3 SmartFlare™ RNA Probe, (Cat. No. SF-471, red), for 16 hours and counter-stained with Hoechst 33342 (blue).

Instead of lysed cells, switch to live cells. And while you're at it, eliminate sample preparation and transfection steps all together. Live cell RNA detection is now possible, in a single incubation step using inert nanoparticle technology to specifically detect native RNA.


And when you're done, the probes exit the cells allowing you to perform downstream analyses with the *same*, unperturbed cells.

**Make the smart change!**

# What's smart

Use a SmartFlare™ probe to analyze gene expression in live cells with no sample prep, cell lysis or toxicity, and you'll see just how smart a tiny particle can be.

- **Live Cell Detection**  
See real-time RNA expression in live cells, giving you more physiologically relevant data.
- **Native RNA**  
No need for amplification, targeting only native, endogenous RNA, avoiding artifacts with traditional amplification methods.
- **No Sample Prep**  
Only a single incubation step, eliminating all sample preparation.
- **Easy Open-platform Detection**  
Detect the fluorescent signals using a platform of your choice: microscope, personal cell analyzer, flow cytometer, cell sorter, etc.
- **No Toxicity, Leaving Cells Unchanged**  
Non-toxic nanoparticles that enter and exit cells via natural mechanisms, allowing the same live cells to be used for downstream assays.

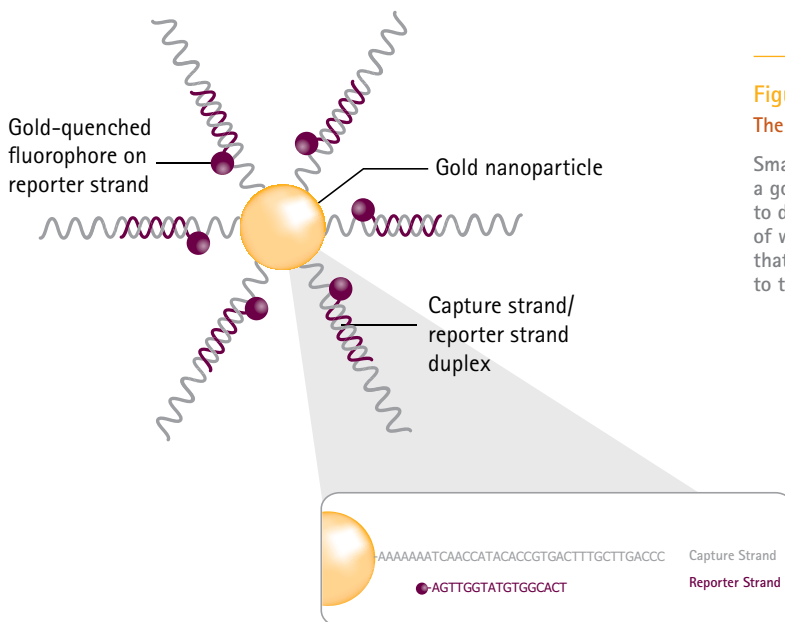

**Figure 1.**  
**The SmartFlare™ Probe Anatomy**

SmartFlare™ probes consist of a gold nanoparticle conjugated to double stranded oligos, one of which contains a fluorophore that is quenched by its proximity to the gold core.

The oligo duplexes are designed with an RNA "capture" sequence and a complementary "reporter" sequence.

# What's simple

Eliminate the complexity and uncertainty associated with other RNA detection techniques. SmartFlare™ RNA Detection Probes make discovery effortless.

- No sample prep
- No RNA amplification
- No transfection reagents
- No cell lysis
- No toxicity


RNA detection often requires laborious sample prep, RNA amplification and detection based on standard curves. SmartFlare™ RNA Detection Probes are endocytosed by live cells using existing cellular endocytosis machinery. Sample prep is unnecessary; simply add SmartFlare™ probes to your cultures, incubate overnight and detect the next day. Over time, the probes will exit the cell, without any adverse effects, allowing for subsequent downstream assays.

## Cell testing protocol

(recommended protocol for testing a cell-probe combination for the first time):

1. Plate cells at desired density (typical example: 30,000 cells in 200  $\mu$ L complete medium in each well of a 96-well plate)
2. Dilute reconstituted SmartFlare™ reagent 1:20 in sterile phosphate-buffered saline to create working solution
3. Add 4  $\mu$ L diluted SmartFlare™ reagent to each well of cells (which are 60-80% confluent)
4. Incubate overnight (16 hours) at appropriate temperature (37 °C), CO<sub>2</sub> (5%) and relative humidity.
5. Detect fluorescence using platform of choice: flow cytometer, imaging cytometer, fluorescence microscope, etc.

# The SmartFlare™ probe's simple mechanism reveals the RNA content of live cells while keeping cells alive and intact.

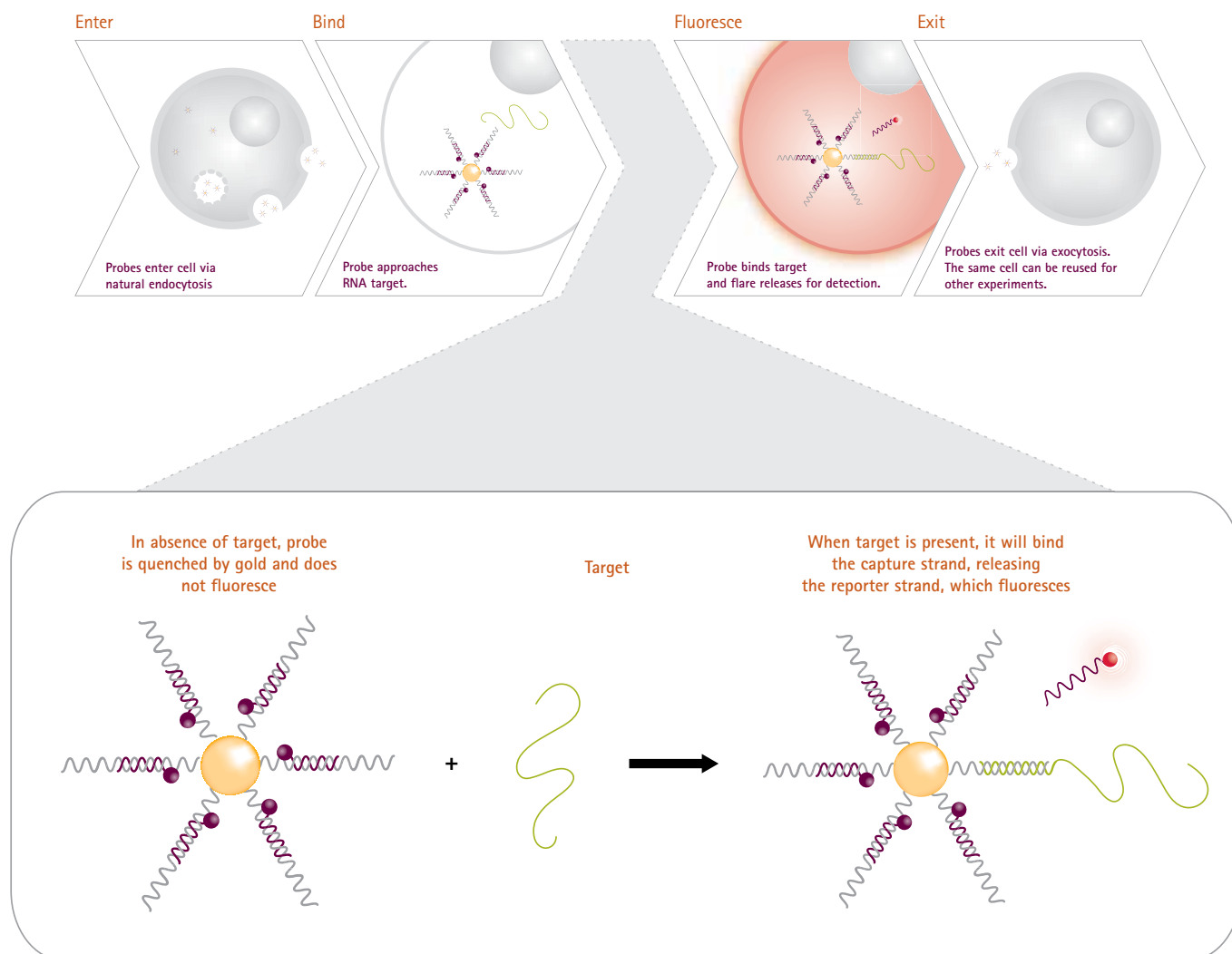

**Figure 2.**

**Enter, bind, fluoresce, exit.** SmartFlare™ probes enter the cell using the cell's own endocytosis process. The probes circulate within the cell and bind to the complementary RNA sequence. This binding event releases a fluorophore, illuminating the cells for detection. Over time, the probe exits the cell, leaving the cell unchanged and free for downstream analyses.

# What's powerful


You can take a peek at your RNA without your cells ever knowing you were there. SmartFlare™ probes enable novel experimental designs that can change how you do science. Imagine the possibilities.

## Applications using SmartFlare™ Probes

- Track changes in RNA dynamically, over time
- Sort cell types that are traditionally difficult to isolate, using intracellular RNA markers
- Quantify miRNAs
- Assess RNA and protein in the same cellular preparations
- Reuse your SmartFlare™ probe-treated cells for downstream assays
- Live cell detection via flow cytometry
- Assess multiple RNA targets simultaneously (multiplex)

## Sort live cells based on RNA content

You can now sort your cells based on levels of RNA targets. After sorting, these cells are alive, unchanged, and ready for further analysis. The ability to detect and separate live cells based on the level of a specific RNA target provides new opportunities to correlate gene expression with cellular functions and to identify rare cell types such as certain tumor cells and cancer stem cells.

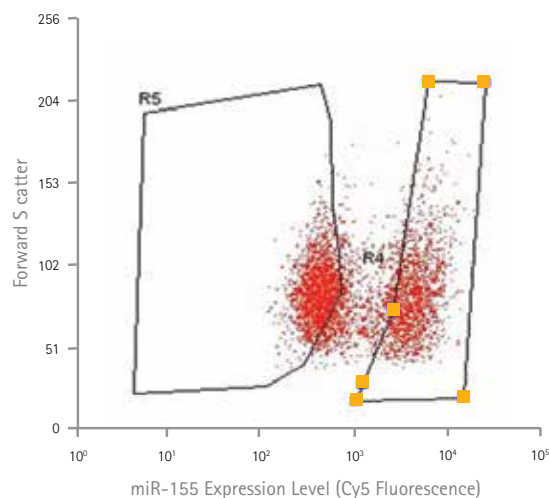


**Figure 3.**  
Sorting based on intracellular RNA markers. Cells were sorted based on miR-155 expression. miR-155 high and miR-155 low cells were isolated and then used for downstream experiments examining protein levels in the same cells.

# What's really happening

Using live cells means you are looking at native, unamplified RNA, which can reveal meaningful links between gene expression profiles and phenotypes. Choose the way you discover what's *really* happening in your cells.

## Visualize RNA using microscopy


Using SmartFlare™ RNA Detection Probes via microscopy is simple. You can see RNA target expression in living cells right under the microscope, track changes in real time, monitor gene expression profiles in mixed cell populations, and get quick, yes/no answers about target expression in living cells.

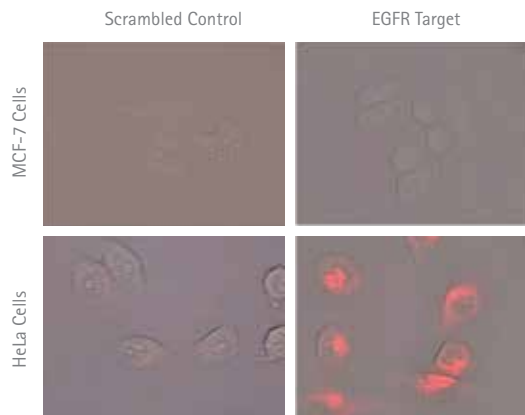

**Figure 4.** Microscopy shows specificity and dynamic range of EGFR SmartFlare™ probe. Varying levels of EGFR expression (red) shown in MCF7 (low expression) and HeLa (high expression) cells when incubated with an EGFR-specific SmartFlare™ probe. No signal is seen in cells treated with a sequence-scrambled control probe (left column).

## Quantify RNA using flow cytometry

Quantifying RNA has never been easier. Because flow cytometry interrogates single cells for fluorescence expression, amplification is not necessary. Get specific, single-cell, relative quantification of RNA in intact cells using EMD Millipore's guava easyCyte™ or Muse™ flow cytometry platforms with SmartFlare™ Probes. Or really enhance your detection capabilities using the Amnis® Imaging Flow Cytometry systems.

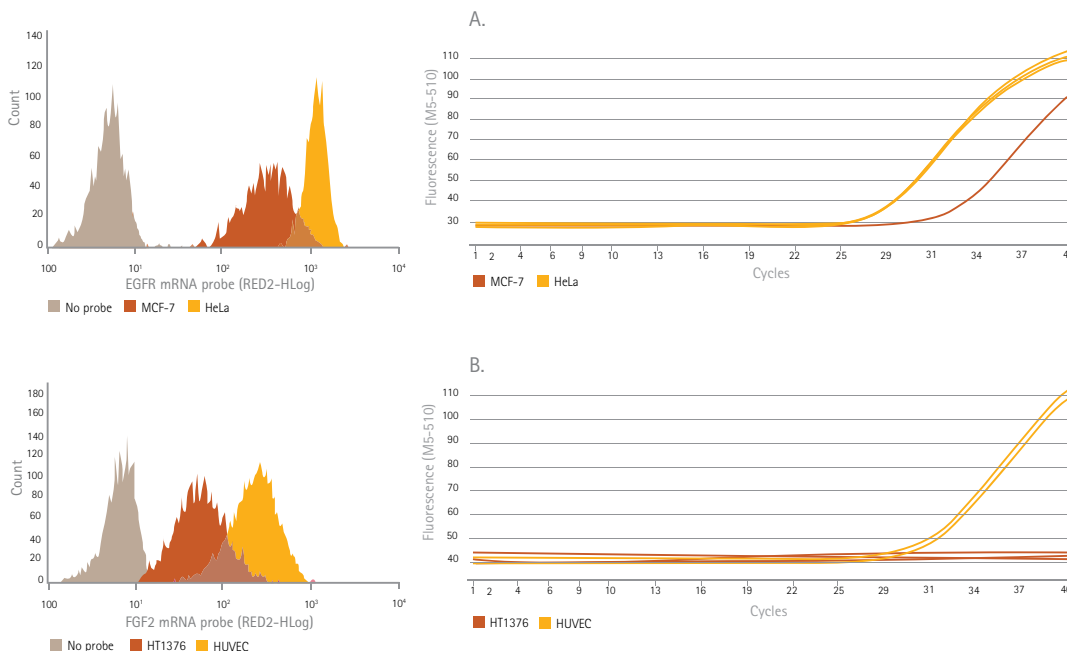

**Figure 5.** Probe detection of mRNA levels correlates to RT-PCR data. Using SmartFlare™ technology to determine the mRNA levels of EGFR mRNA (A) in HeLa and MCF-7 cells as well as FGF2 mRNA (B) in HUVEC and HT1376 cells. RNA levels as detected by flow cytometry show good correlation to trends in cycle numbers and RNA levels determined by quantitative RT-PCR.

# What's the best detection platform

SmartFlare™ RNA Detection Probes are compatible with almost any fluorescence detection system. After a simple overnight incubation, you are ready to analyze your cell's RNA content on the platform of your choice.

Choose from a wide range of fluorescence detection platforms based on your research needs.

	Microscopy	Personalized Cell Analyzer	Flow Cytometry	Imaging Flow Cytometry	Cell Sorter
Live Cell Analysis	✓	✓	✓	✓	✓
Visualization	✓			✓	
Relative Quantitation		✓	✓	✓	✓
Single Cell Analysis	✓	✓	✓	✓	✓
Sorting Capabilities					✓
96-well Format			✓	✓	✓




# EMD Millipore offers powerful and innovative options for single cell fluorescence detection.

## Muse™ Cell Analyzer

- Compact and affordable
- Intuitive software and touchscreen interface allow for easy operation - no flow expertise necessary
- Green laser for flow-based detection of Cy3 fluorophore
- Optimized assays for cell health, cell signalling and immunology

[www.millipore.com/muse](http://www.millipore.com/muse)

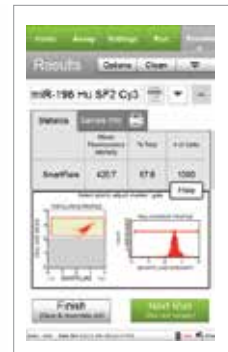


Simplified analysis of miR-196 in single cells with SmartFlare™ probes and the Muse™ Cell Analyzer.

## guava easyCyte™ Benchtop Flow Cytometer

- Automated sampling for 96-well plates
- No sheath fluid—microcapillary system lets you use less sample and less reagents for analysis
- Powerful InCyte™ software with capabilities for heatmapping and real-time adjustment of data analysis

[www.millipore.com/guava](http://www.millipore.com/guava)

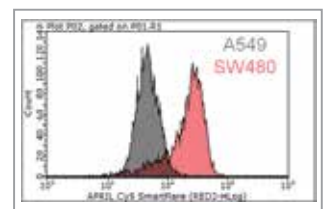


Population distribution with respect to APRIL expression in two cell lines as determined using the guava easyCyte™ Benchtop Flow Cytometer and SmartFlare™ probes.

## Amnis® Brand Imaging Flow Cytometers

- Combines the power of digital fluorescence microscopy with the speed and sensitivity of flow cytometry
- Perform phenotypic and functional studies at the same time using up to seven lasers and 12 images per cell

[www.amnis.com](http://www.amnis.com)


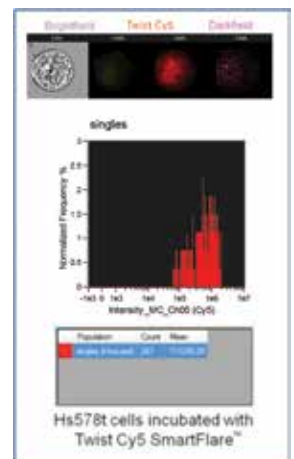


Imaging flow cytometry using the Amnis® brand instrument and SmartFlare™ probes reveals the pattern and intensity of Twist expression in individual Hs578t cells.



# What's the first step

Interested in testing SmartFlare™ probes in your system? Start with our controls. We recommend starting with an Uptake and Housekeeping control to see how SmartFlare™ RNA detection probes can help your research. Whether you're an experienced user or just interested in performing some initial proof of concept experiments, we have the right control for you. Be confident in your results by having the proper SmartFlare™ controls in your experiments.

## A. Uptake Control

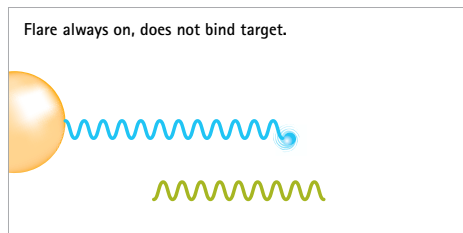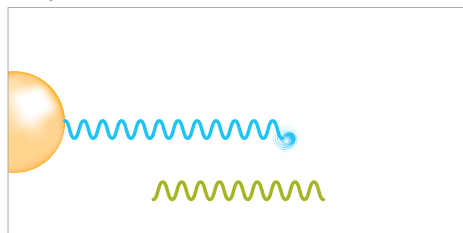

## B. Scramble Control

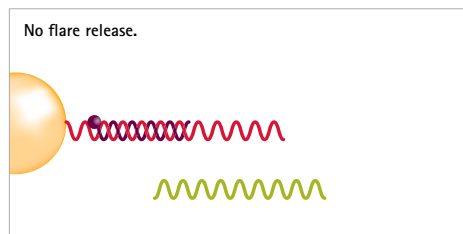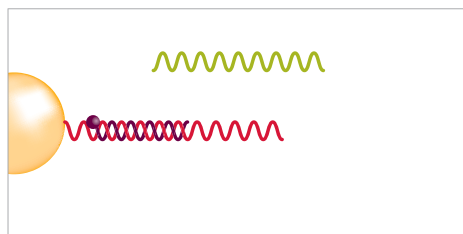

## C. Housekeeping Control

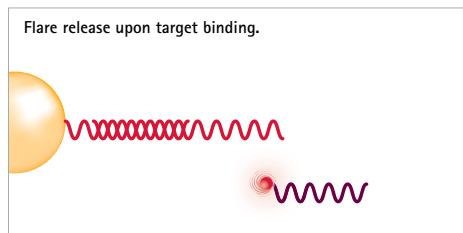
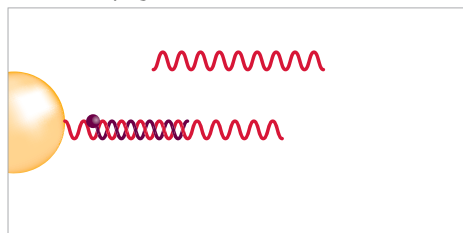


Figure 6.

### SmartFlare™ Control:

- A. Uptake Control – These controls are always fluorescing and are a good way of confirming the presence and relative amount of SmartFlare™ probes that enter your cell of interest.
- B. Scramble Control – These controls help you account for any background fluorescence within the cell.
- C. Housekeeping Control – These controls target commonly-expressed genes, are great for introductory experiments and serve as positive controls for your target SmartFlare™ probes.

## Related Products - SmartFlare™ Controls

Description	Qty.	Cat. No.
<b>Housekeeping Controls</b>		
18S Hu, Ms, Rt-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-143
18S Hu, Ms, Rt-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-142
Beta Actin Hu-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-145
Beta Actin Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-144
Cyclophilin A Hu-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-150
Cyclophilin A Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-139
Cyclophilin A Ms-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-129
Cyclophilin A Ms-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-128
GAPDH Hu-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-126
GAPDH Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-136
GAPDH Ms-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-125
GAPDH Ms-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-138
GAPDH Rt-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-141
GAPDH Rt-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-140
HPRt1 Hu-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-133
HPRt1 Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-132
HPRt1 Ms-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-135
HPRt1 Ms-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-134
HPRt1 Rt-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-112
HPRt1 Rt-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-189
RNA Pol II Hu-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-131
RNA Pol II Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-130
RNA Pol II Ms-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-123
RNA Pol II Ms-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-124
RPLP0 Hu-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-482
RPLP0 Ms-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-187
RPLP0 Ms-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-199

Description	Qty.	Cat. No.
<b>Uptake Controls</b>		
Uptake-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-114
Uptake-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-137
<b>Scramble Controls</b>		
Scramble-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-103
Scramble-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-102
miRNA Scramble-Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-147
miRNA Scramble-Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-146
<b>Uptake/Scramble Combination Controls</b>		
Scramble Cy3+Uptake Cy5 SmartFlare™ RNA Detection Probe	250 reactions	SF-105
Scramble Cy5+Uptake Cy3 SmartFlare™ RNA Detection Probe	250 reactions	SF-104

Please visit our website for the complete list of SmartFlare™ RNA detection probes. More performance data and application notes can also be found on our website.

[www.millipore.com/smartflare](http://www.millipore.com/smartflare)

To learn more and order online, visit:

[www.millipore.com/smartflare](http://www.millipore.com/smartflare)



## To Place an Order or Receive Technical Assistance

In the U.S. and Canada, call toll-free 1-800-221-1975

For other countries across Europe and the world, please visit: [www.emdmillipore.com/offices](http://www.emdmillipore.com/offices)

For Technical Service, please visit: [www.emdmillipore.com/techservice](http://www.emdmillipore.com/techservice)

## Get Connected!

Join EMD Millipore Bioscience on your favorite social media outlet for the latest updates, news, products, innovations, and contests!

 [facebook.com/EMDMilliporeBioscience](https://facebook.com/EMDMilliporeBioscience)

 [twitter.com/EMDMilliporeBio](https://twitter.com/EMDMilliporeBio)



[www.emdmillipore.com/offices](http://www.emdmillipore.com/offices)