

Technical Bulletin

Oxaloacetate Assay Kit

Catalogue number **MAK515**

Product Description

Oxaloacetate (OAA) is an intermediate in the citric acid cycle and participates in gluconeogenesis. OAA is formed by the oxidation of malate, by deamidation of aspartate or by condensation of CO₂ with pyruvate or phosphoenolpyruvate.

The Oxaloacetate Assay Kit provides a simple, direct and automation-ready procedure for measuring oxaloacetate concentration. OAA is converted into pyruvate which is then oxidized with the conversion of the dye into a colored and fluorescent form. The color intensity of the oxidized dye at 570 nm or fluorescence intensity at $\lambda_{\text{ex}} = 530 \text{ nm}$ / $\lambda_{\text{em}} = 585 \text{ nm}$ is directly proportional to the oxaloacetate concentration in the sample.

The linear detection range of the kit is 7 to 400 μM oxaloacetate for the colorimetric assays and 1 to 40 μM for the fluorometric assays. The kit is suitable for oxaloacetate determination in plasma, serum, tissue and culture media.

Components

The kit is sufficient for 100 colorimetric OR Fluorometric assays in 96-well plates.

- Developer 10 mL
Catalogue Number MAK515A
- ODC Enzyme 120 μL
Catalogue Number MAK515B
- Dye Reagent 120 μL
Catalogue Number MAK515C
- Oxaloacetate Standard 1 vial
Catalogue Number MAK515D

Equipment Required but Not Provided

- Pipetting devices and accessories (for example, multichannel pipettor)
- Multiwell plate reader
- Clear flat-bottom 96-well plates for colorimetric assay or black flat-bottom 96-well plates for fluorometric assay. Cell culture or tissue culture treated plates are not recommended.
- microcentrifuge tubes
- Dounce tissue grinder set (Catalogue Number D9063 or equivalent)
- Refrigerated microcentrifuge capable of $\geq 14,000 \text{ rpm}$.

Precautions and Disclaimer

For R&D use only. Not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

The kit is shipped on wet ice. Store components at $-20 \text{ }^\circ\text{C}$.

Preparation Instructions

Briefly centrifuge small vials prior to opening. Equilibrate all components to room temperature prior to use.

Colorimetric Procedure

All Samples and Standards should be run in duplicate.

Sample Preparation

Tissue and Cell Samples

1. Homogenize the tissue and cell Samples (2×10^6) in 100 μL PBS.
2. Centrifuge at 14,000 rpm for 5 minutes.
3. Deproteinize the supernatant using a 10 kDa spin filter (Catalogue UFC8010 or equivalent).

Serum and Plasma Samples

Samples must be deproteinized and an internal Standard should be used.

Cell Culture Media

Avoid media with high pyruvate concentrations when assaying Cell Culture Media Samples.

(For example, DMEM, L-15, F12, etc.)

Note: If using an internal Standard, Samples will need three separate reactions.

Sample Blank

1. Transfer 20 μL of each Sample to separate wells.
2. To each Sample Blank well, add 5 μL of purified water.

Internal Standard

(Required for serum and plasma Samples)

1. Transfer 20 μL of each Sample to separate wells.
2. Prepare 500 μL 80 μM OAA Standard by mixing 100 μL 400 μM OAA standard and 400 μL purified water.
3. For colorimetric and fluorometric assays, to each Internal Standard well, add 5 μL of 80 μM OAA Standard.

Add 20 μL of each Sample and 5 μL of dH_2O to two separate wells in a 96 well plate.

Standard Curve Preparation

1. Prepare 10 mM stock by dissolving it with 100 μL purified water.
2. Prepare a 400 μM OAA Standard by diluting 20 μL of the 10 mM Standard with 480 μL purified water.
3. Further, dilute Standards in 1.5mL centrifuge tubes as mentioned in Table 1.

Note:

1. Keep standard cold and store at $-20\text{ }^\circ\text{C}$.
2. Reconstituted OAA Standard should be used within 2 weeks.
3. If assaying culture media with phenol red, dilute the Oxaloacetate Standard in culture media.

Table 1.

Dilution of Standard

No.	400 μM OAA Standard	Purified H_2O	Oxaloacetate (μM)
1	100 μL	0 μL	400
2	60 μL	40 μL	240
3	30 μL	70 μL	120
4	0 μL	100 μL	0

4. Transfer 20 μL Standards into separate wells of a clear flat-bottom 96-well plate.

Working Reagent Preparation

1. Mix enough reagents for the number of assays to be performed. For each Standard and Sample well, prepare the Working Reagent according to Table 2.

Table 2.

Preparation of Working Reagents

Reagent	Working Reagent	Blank Working Reagent
Developer	85 μL	85 μL
ODC Enzyme	1 μL	-
Dye Reagent	1 μL	1 μL

2. Transfer 80 μL Working Reagent into each reaction well. Add 80 μL Blank Working Reagent to the Sample Blank Wells. Tap plate to mix.

Measurement

1. Incubate the plate, protected from light, for 15 minutes at room temperature.
2. Read the optical density at 570 nm.

Fluorometric Procedure

1. Dilute the Standards prepared in Colorimetric Procedure 1:10 in purified water. If an internal Standard is used, use the same concentration as described in the Colorimetric Procedure (for example, 5 μL of 80 μM OAA).
2. Transfer 20 μL of the Standards and 20 μL of Sample (2 wells per Sample if a Standard curve is used; 3 wells per sample if an internal Standard is used, see Colorimetric Procedure) into separate wells of a black 96-well plate.
3. Add 80 μL of appropriate Working Reagent (see Colorimetric Procedure) to each well. Tap the wellplate to mix.

Measurement

1. Incubate for 15 minutes at room temperature.
2. Read fluorescence at $\lambda_{\text{ex}} = 530 \text{ nm}$ and $\lambda_{\text{em}} = 585 \text{ nm}$.

Results

1. Calculate ΔOD or ΔF by subtracting the reading (OD or fluorescence intensity F) of Standard #4 (Blank) from the remaining Standard reading values.
2. Plot the ΔOD or ΔF against Standard concentrations and determine the slope of the Standard curve.
3. Calculate the Oxaloacetate Concentration of samples using the below given equation:

$$\text{Oxaloacetate } (\mu\text{M}) = \frac{R_{\text{Sample}} - R_{\text{Blank}}}{\text{Slope}(\mu\text{M}^{-1})} \times \text{DF}$$

If an internal standard was used, the sample Oxaloacetate concentration is computed as follows:

$$\text{Oxaloacetate } (\mu\text{M}) = \frac{R_{\text{Sample}} - R_{\text{Blank}}}{R_{\text{Standard}} - R_{\text{Sample}}} \times \frac{80}{4} \times \text{DF}$$

Where:

- R_{Sample} = OD or fluorescence intensity (F) reading of Sample
- R_{Blank} = OD or fluorescence intensity (F) reading of Blank
- 80 = Concentration in μM of Internal Standard
- 4 = The volume of the Internal Standard is 4 \times lower than the Sample volume. In order to correct the calculation, the Internal Standard concentration is divided by 4.

DF = Sample dilution factor (DF = 1 for undiluted samples)

Conversions: μM oxaloacetate equals 13.1 mg/L, 0.00131% or 13.1 ppm.

Figure 1.

Typical Colorimetric Oxaloacetate Standard Curve.


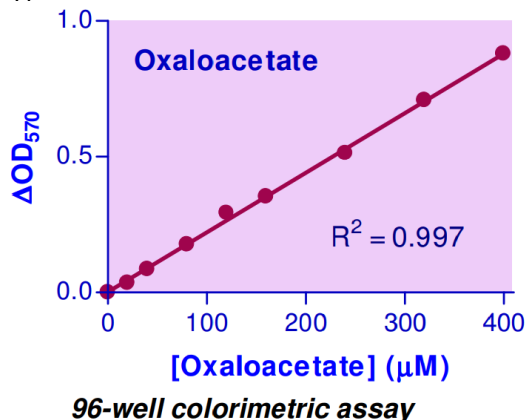


Figure 2.

Typical Fluorometric Oxaloacetate Standard Curve

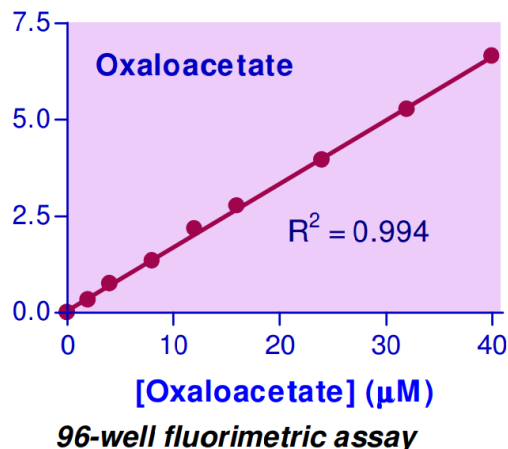

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the U.S. and Canada.

MilliporeSigma, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2024 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

mak515pis Rev 08/24

**MILLIPORE
SIGMA**